



 **CHATTANOOGA GROUP**
Moving Rehabilitation Forward™

TX Traction Products



New innovation for a classic therapy

TX® Traction System

4759 – TX Traction Unit

Standard Accessories


48082 – Patient Interrupt Switch

27465 – Patient Data Cards (5x)


USER PROTOCOLS



CLINICAL LIBRARY

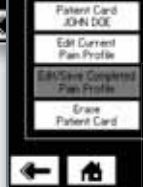

PAIN MAP

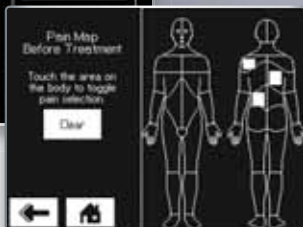

PAIN MAP



In 1979 Chattanooga Group introduced the TX line of clinical traction devices to the physical medicine market. Today, Chattanooga Group is the world leader in clinical traction devices, tables and accessories. Using the most advanced technology, the new TX traction unit is years ahead of any traction unit currently available on the market today. With innovative features and an easy to use digital touch screen interface that makes set up simple, the redesigned TX traction unit is the first true innovation to this classic therapy in over 25 years.

TX Innovative Features


- Easy to operate
- Digital Touch Screen Interface
- 270° Pivoting Interface
- Patient Pain Mapping
- Patient Data Card records pre- and post-treatment pain profiles, treatment session information and patient information
- 10 user-defined protocols
- Real-time digital treatment tracking
- User-defined Hold, Rest and Treatment times
- Intermittent, Static and Cyclic traction programs
- Progressive and Regressive settings
- Traction tension parameters 0 – 200 lbs.
 - » Patient Interrupt Switch
 - » Cervical maximum poundage warning
 - » Audible signal at end of treatment or when Patient Interrupt Switch is activated

Technical Specifications

Mains Power: 100 – 240 V, 50/60 Hz
 Power Consumption: 75 W
 Current Consumption: 3.2 Amps Max
 Weight: 30 lbs (14 kg)
 Dimensions: 17.5" L x 9.5" W x 17.5" H
 (45 cm x 24 cm x 45 cm)
 Electrical Safety Class: Class 1, Type B
 Safety Tests: IEC/UL/EN: 60661-1, 60601-1-2;
 CAN C22.2 No. 601.1-M90 w/AZ;
 cETLus Listed



TRACTION TENSION PARAMETERS	MIN - MAX	INCREMENTS
TRACTION PERIOD	1 - 99 MINUTES	1 MINUTE
HOLD PERIOD	0 - 99 SECONDS	1 SECOND
REST PERIOD	0 - 99 SECONDS	1 SECOND
TRACTION TENSION	0 - 200 LB (0 kg/0 N) - (90kg/890 N)	1 LB (1 kg/5 N)
PROGRESSIVE STEPS	1 STEP - 9 STEPS	1 STEP
REGRESSIVE STEPS	1 STEP - 9 STEPS	1 STEP

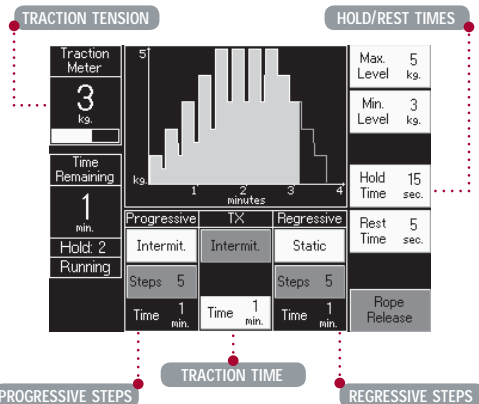

Patient Data Management System

2790 – Patient Data Management System (PDMS)
 Includes: software, card reader and 5 patient data cards

27465 – Patient Data Cards (25/pack)

Mobile Traction Stand

3600 – Mobile Traction Stand

Patient Data Management System (PDMS)



PDMS is a unique optional software program that allows the clinician to organize, store, display and print patient traction therapy information from the TX device to the PC. This system which includes software and a card reader for your PC can be used with the patient data cards (5 included with PDMS) to document treatments.

Mobile Traction Stand

- Fully adjustable
- Allows the TX Traction device to be used with most tables
- Traction unit sold separately





ADP-400 Traction Table

- Four sections
- Friction free gliding lumbar section
- Turret mounted traction pedestal
- Head and foot sections incline up to 90°
- Head section lowers to 25°
- 400 lb (181 kg) lifting capacity
- Adjustable height range between 21.5" to 41.5" (55 cm to 105 cm)
- Standard hand control
- Optional foot control
- Optional retractable caster
- Available in all upholstery colors*
- Traction unit sold separately



Technical Specifications

Mains Power: 120 V, 50 Hz
 Current: 1.6 A
 Power Input: 200 W
 Weight: 140 lbs (64 kg)
 Shipping Weight: 200 lbs (91 kg)
 Dimensions: 77"L x 27"W x 21.5" to 41.5"H
 (196 cm x 69 cm x 55 cm to 105 cm)
 Electrical Class: Class 1, Type B

ADP-400 Traction Table

ADP400 – ADP-400 Traction Table

TTET-200 Traction Table

6870 – TTET-200 Traction Table

TTFT-200 Traction Table

6840 – TTFT-200 Traction Table

Standard Colors



TTET-200 Traction Table

- Two sections
- Friction free gliding lumbar section
- Traction unit pedestal
- 400 lb (181 kg) lifting capacity
- Adjustable height range between 19" to 38" (48 cm to 96 cm)
- Standard hand control
- Optional foot control
- Available in all upholstery colors*
- Traction unit sold separately



Technical Specifications

Mains Power: 120 V, 50 Hz
 Current: 1.6 A
 Power Input: 200 W
 Weight: 140 lbs (64 kg)
 Shipping Weight: 200 lbs (91 kg)
 Dimensions: 75"L x 28"W x 19" to 38"H
 (191 cm x 71 cm x 48 cm to 96 cm)
 Electrical Class: Class 1, Type B

Designer Colors*



RhinoHyde® Colors*



TTFT-200 Traction Table

- Two sections
- Traction unit pedestal
- 350 lb (159 kg) working load capacity
- Fixed height of 32" (81 cm)
- Available in all upholstery colors*



*Additional charges apply for Designer and RhinoHyde colors.

This color chart is produced within the limitations of lithographic color printing and is intended as a general guide only. RhinoHyde® Products are a registered trademark of Tandem Products, Inc.



Saunders Cervical Traction

7040 – Saunders Cervical Traction Device: includes clevis for attachment to TX traction device (48039)

Carpal Trac

1450 – Carpal Trac: includes Clevis for attachment to TX traction device (48039)

12450 – Replacement Wrist Strap

TX Clevis

48039 – TX Clevis

TXS-1 Flexion Stool

7050 – TXS-1 Flexion Stool: adjustable height from 13-20" (33 cm - 51 cm)

Overdoor Traction Kit

4390 – Overdoor Traction Kit: includes overdoor pulley, water bag, and halter

TXA-1 Accessory Package

1440 – TXA-1 Accessory Package Includes:

- One adjustable cervical traction halter (1405)
- One heavy duty pelvic traction set (1433)
- Two TX pillows (1301) and covers (1341)
- One 17" (43 cm) spreader bar (60030)

Heavy Duty Pelvic Traction

1433 – Heavy Duty Pelvic Traction Set (Universal fit to 52") Includes:

- Heavy Duty Thoracic Restraint (1435)
- Heavy Duty Pelvic Belt (1434)
- Thoracic Restraint Straps (1432)

1435 – Heavy Duty Thoracic Restraint

1434 – Heavy Duty Pelvic Belt

1432 – Thoracic Restraint Straps

11261 – Replacement Strap with Buckle (not shown). For use with heavy duty thoracic restraint and heavy duty pelvic belt

TX Pillow

1301 – TX Pillow (without cover)

TX Pillow Cover

1341 – TX Pillow cover

Double Pull Pelvic Belt (not shown)

1412 – Medium

1413 – Large

1414 – X-Large

TX® Traction Harnesses

Cervical Traction Halter (A)

1403 – Medium

1404 – Large

1405 – Adjustable

Universal Cervical Traction Halter (B)

1401 – Universal Cervical Traction Halter

Universal Thoracic Restraint (C)

1431 – Universal Thoracic Restraint

Universal Single Pull Pelvic Belt (D)

1421 – Universal Single Pull Pelvic Belt

Naugahyde® is a registered trademark of Uniroyal Engineered Products, LLC.

Chattanooga Group Exclusive

Saunders Cervical Traction

This unit pulls at base of the occiput for more effective cervical traction. It adapts easily to all Chattanooga Group traction devices. One size fits most patients.



Note:

Clevis included to mount to traction device.



TXS-1 Flexion Stool

Adjustable height feature allows therapist to place the patient's lumbar spine into neutral or flexed position for increased comfort. Upholstery can be color-coordinated to match traction table.

Carpal-Trac

This traction accessory is designed to alleviate pressure on the soft tissue structures of the wrists. A comfortable, inter-locking wrist strap allows therapeutic traction without compressing carpal tunnel soft tissue.



Overdoor Traction Kit

This effective, over-the-door cervical traction kit contains everything needed for patient-administered flexion or extension cervical traction.



TXA-1 Traction Accessory Package

For pelvic or lumbar traction applications up to 200 lbs (91 kg). Includes pelvic belt, thoracic restraints and thoracic straps with durable and comfortable Naugahyde® pads to reduce slippage.



Traction Harnesses

TX® Cervical Halters, Thoracic Restraints, and Pelvic Belts: The thoracic restraint and pelvic belt accommodate pelvic or lumbar applications up to 200 lbs (91 kg).

